





PIXIS: 2K

The PIXIS series from Princeton Instruments (PI) are fully integrated, low noise cameras designed for quantitative scientific imaging optical spectroscopy applications. Designed utilizing PI's exclusive XP cooling technology, PIXIS offer thermoelectric cooling down to -75°C, while the all-metal, hermetically sealed design, with the industry's only lifetime vacuum guarantee, ensures maintenance free operation. Choose either back-illuminated (B/BUV) or eXcelon versions of the 2048 x 512 CCD, for optimized performance from UV to NIR. High QE and ultra low-noise electronics make the PIXIS: 2K series of cameras ideal for demanding, low light applications such as Raman spectroscopy, LIBS, absorbance, emission, fluorescence and reflectance spectroscopy. Princeton Instruments' eXcelon technology delivers the highest sensitivity in the NIR while suppressing the etaloning that occurs in standard back illuminated CCDs. Dual speed operation at 100kHz or 2MHz allow these cameras to be used for steady state as well as fast kinetics studies.

FEATURE	BENEFITS	
eXcelon® technology	Increases detector sensitivity while suppressing etalon interference fringes observed in the NIR with conventional back-illuminated devices.	
Permanent vacuum	Our all-metal brazed seals eliminate the outgassing that occurs with epoxy seals and allow us to offer the industry's only lifetime vacuum guarantee.	
Deep thermoelectric cooling	Low dark noise allows detection of faint signals; CoolCUBE II, a compact room temperatur coolant circulator, is available for vibration sensitive environments.	
Single fused silica vacuum window	Minimizes reflection losses from the UV to the NIR; Optional AR coating and wedge windows are available.	
Optional UV phosphor coating	Enhances sensitivity throughout the UV to below 200 nm.	
2048 x 512 imaging array, 13.5μm x 13.5μm pixels	CCD format provides 6.9 mm height for rapid spectral acquisition; 13.5 µm pixel size offers superior resolution to industry-standard detectors with a 1024 pixel format.	
100 kHz and 2 MHz digitization rates	Dual-speed digitization allows complete freedom to select between slow operation for low noise and highest SNR or fast operation for rapid image acquisition.	
Dual amplifiers with software-selectable system gains	ins High sensitivity amplifier reduces noise for weak signals; High capacity amplifier increases dynamic range.	
USB2.0 data interface	Plug-and-play operation with desktops or laptops; Optional fiber optic interface for remote operation; Easy OEM integration.	
Optional: LightField® (for Windows 8/7, 64-bit) Or WinView/Spec (for Windows 8/7/XP, 32-bit)		
Compatible with Windows 8/7/XP, and Linux; Oftware development kits (SDKs) Compatible with Windows 8/7/XP, and Linux; Universal programming interfaces for easy custom programming.		

Raman Spectroscopy, LIBS, Absorbance, Emission, Fluorescence and Reflectance Spectroscopy



SPECIFICATIONS



	enabled			
	PIXIS: 2KB_eXcelon	PIXIS: 2KB/BUV		
Features	Back-illuminated CCD with eXcelon technology. Enhanced sensitivity in the UV and the NIR with low etaloning.	Back-illuminated CCD. High sensitivity in both the visible and UV regions. Special BUV version offers the highest sensitivity in the UV region.		
Dark current @ -75°C (e-/p/sec)	0.001 (typical)	0.001 (typical)		
CCD format	2048 x 512 imaging pixels; 13.5 x 13.5-μm pixels with 100% fill factor			
Imaging area	27.6 x 6.9 mm (optically centered)			
Optical mount	Princeton Instruments Acton spectrometer adapter with optional shutter			
Deepest cooling temperature	-75°C guaranteed			
Thermostating precision	±0.05°C			
Cooling method	Thermoelectric air or liquid cooling (CoolCUBE II required)			
Spectrometric well capacity: High Sensitivity High Capacity	250 ke- (typical), 150 ke- (min) 750 ke- (typical), 600 ke- (min)			
ADC speed/bits	100kHz/16-bit and 2MHz/16-bit			
System read noise @100 kHz @2 MHz	3.5 e- rms (typical), 6 e- rms (max) 14 e- rms (typical), 20 e- rms (max)			
Vertical shift speed	15.2 μsec/row			
Non-linearity	<1% @ 100 kHz			
Software selectable gains	1.5, 3, 6 e- (high sensitivity); 3, 6, 12 e- (high capacity); available at all speeds			
Operating systems supported	Windows 8/7/XP (32-bit), Windows 8/7 (64-bit) and Linux			
Data interface	USB2.0			
I/O signals	Two MCX to BNC connectors for programmable frame readout, shutter, trigger in			
Operating environment	+5 to +30°C non-condensing			
Certification	CE			
Dimensions / Weight	16.59 cm (6.53") x 11.81 cm (4.65") x 11.38 cm (4.48") (L x W x H) / 2.27kg (5lbs)			

All specifications are subject to change.

SPECTRAL RATE

@ 100 kHz Full Vertical Binning (FVB)		35 spectra/sec
	@ 2 MHz Full Vertical Binning (FVB)	90 spectra/sec

QE DATA



NOTE:

Wavelength (nm)


Graph shows typical Quantum Efficiency (QE) data measured at $+25^{\circ}$ C. QE decreases at normal operating temperatures. For the best results for your application, please discuss the specific parameters of your experiment with your sales representative.

VACUUM WINDOW AR COATINGS

NOTES:

- Standard anti-reflection (AR) coating options shown on graph
- Designed by Acton Optics, our BBAR coating offers unmatched performance for 400 nm - 1100 nm
- Custom wedge window options and other AR coatings are also available


Contact your local sales representative for more information

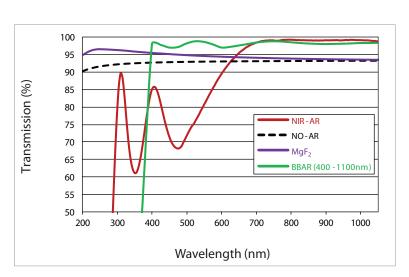




eXcelon Performance

QE Improvement (B_eXcelon vs. B)

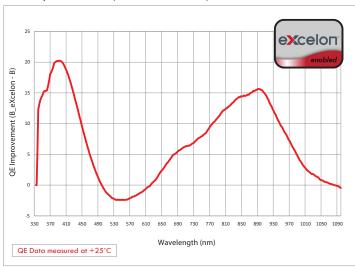


B_eXcelon provides superior QE over the standard back illuminated ("B") version in the UV-NIR range.

Etalon Oscillations (B_eXcelon vs. B)

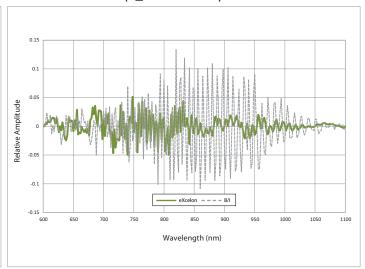


B_eXcelon provides significantly lower etaloning (unwanted fringes) compared to standard back illuminated ("B") version.

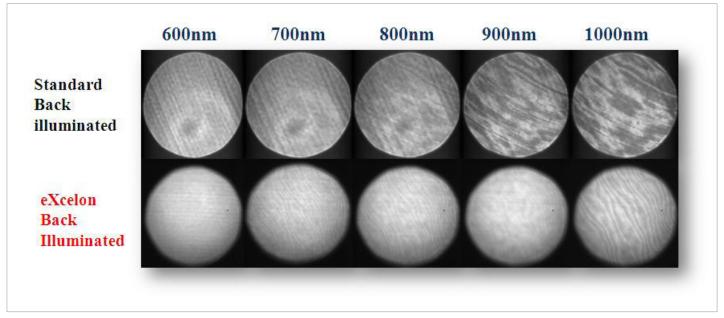


Data taken with white light source through a monochromator, comparing etaloning performance of eXcelon vs. back-illuminated CCDs.

PIXIS DRAWINGS (AIR COOLED)

PIXIS with Spectroscopy Mount

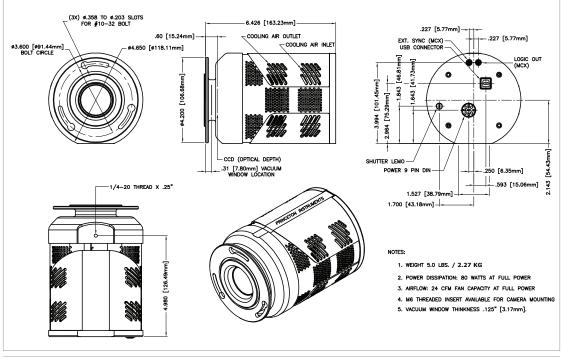

PIXIS with Spectroscopy Mount; with Shutter

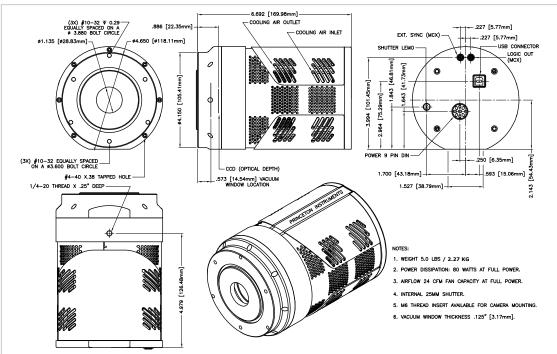

PIXIS DRAWINGS (LIQUID COOLED)

